EN

DE

TR

FR

ES

IT

RO

PL

RU

CRO

BIH

GE

GR

PT

AR

**Zilan®**

*Digital Body Scale*

ZLN8979

ZLN8980

**INSTRUCTION MANUAL**

[www.zilan.com.tr](http://www.zilan.com.tr)

**2 X AAA 1.5 V  
Battery**

PLEASE NOTE THAT THE DEVICE MANUAL CAN BE CREATED FOR MORE THAN ONE MODEL.  
INFORMATION YOUR DEVICE CAN NOT BE CONTAINING THE FUNCTION LISTED IN THE MANUAL.

## GENERAL SAFETY INSTRUCTIONS & WARNINGS

- If the surface of the scale becomes dirty, wipe it with a clean cloth and alcohol or glass cleaner. Keep it away from soap, water, abrasive liquids, high heat, or extreme cold.
- Check the battery if the scale is not functioning properly, and replace it if necessary.
- Store the batteries in a place out of reach of young children and take precautions to prevent accidental swallowing. Seek medical help immediately if a battery is swallowed.
- If the scale will not be used for a long period, remove the battery and store it separately.
- If the batteries are damaged or leaking, dispose of them safely and avoid contact.

## FEATURES

- Maximum weight capacity: 3 kg - 180 kg
- High-precision measurement sensors (Division: 1D = 50 g)
- Unit conversion option (kg/lb/st)
- Auto power-off function
- Low battery and overload indicators

## INSTALLING THE BATTERIES

- Open the battery compartment cover.
- Insert 2 AAA 1.5V batteries, ensuring correct polarity, and close the cover.
- When the battery level is low, a warning symbol will appear on the display. Replace the batteries immediately.

## INSTRUCTIONS FOR USE

- Place the scale on a flat surface. (Measurements taken on uneven surfaces or soft surfaces like carpets and rugs may be inaccurate.)
- Press the small button located at the bottom of the scale to select the weight unit (kg, lb, or st).
- Step onto the scale barefoot and balance your weight.
- The measurement will appear on the display approximately 3 seconds after stepping on the scale.
- The scale will automatically turn off after 8 seconds.

### Notes

- Ensure the scale resets to zero before each use.
- Use the scale only for personal weight measurement.
- Remove the batteries during long-term storage.

## MAINTENANCE AND CLEANING

- Wipe the scale with a damp cloth. Do not use detergents or abrasive cleaners.
- Store the device in an upright position to protect the sensitive sensors.

## OKOK SETUP

### S1 - Install App

- Switch on Bluetooth
  - Phone or Tablet connects to Net (Wi-Fi, 3G or 4G)
- Use your device to scan QR Code or search OKOK International from App Store/Google Play to install App.


### S2 - Pair the Devices, Register and Complete User Information

Open APP, begin to weigh and pair.

Register as a user.

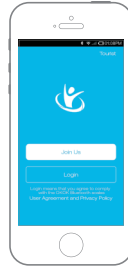Complete personal information.

### S3 - Weighing

Enter into "Dynamic Weight" interface.

Stand on the scale and synchronize weight data onto APP.


### S4 - View the Test Result

Enter into "Weight Trend" interface.

View test report.
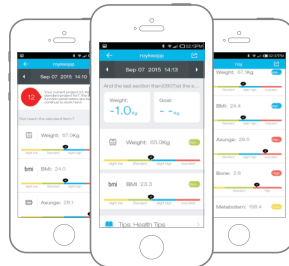
### Correct Disposal of this product

- This marking indicates that this product should not be disposed with other household wastes throughout the EU.
- To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.
- To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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