FLORIA®

Make it Easy

Instruction Manual

Mechanical Body Scale

ZLN1693

EN

DE

TR

FR

ES

IT

RO

PL

RU

CRO

BIH

GE

AL

BG

GR


HU

AR

PLEASE NOTE THAT THE DEVICE MANUAL CAN BE CREATED FOR MORE THAN ONE MODEL.
INFORMATION YOUR DEVICE CAN NOT BE CONTAINING THE FUNCTION LISTED IN THE MANUAL.

MECHANICAL BODY SCALE

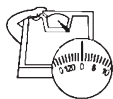
Adjusting to "0"

To adjust the dial to zero, turn the taring wheel at the base of the scale until the zero aligns with the pointer.

Correct Placement

Place the scale on a solid, level floor, as a surface such as a carpet can seriously affect the reading.

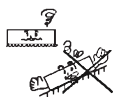Correct Usage

Stand squarely and gently on the scale without leaning forward or backward, distributing your weight evenly between your feet. Do not jump on the scale.


Maintenance

Use a soft cloth slightly moistened with a mild detergent solution. Do not use strong solvents such as benzene or thinner. Oiling the internal mechanism is not required.

Scope of Application

This scale is designed for measuring body weight and conforms to international standards for similar products. It is not intended for commercial weighing purposes.

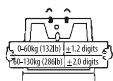Accuracy

If the scale is used properly according to the above instructions, its accuracy will be within ± 1.2 digits for body weights from 0 to 60 kilograms and within ± 2.0 digits for body weights from 60 to 130 kilograms.

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