FLORIA®

Make it Easy

Instruction Manual

Digital Body Scale

ZLN8962

EN

DE

TR

FR

ES

IT

RO

PL

RU

CRO

BIH

GE

GR

PT

AR

PLEASE NOTE THAT THE DEVICE MANUAL CAN BE CREATED FOR MORE THAN ONE MODEL.
INFORMATION YOUR DEVICE CAN NOT BE CONTAINING THE FUNCTION LISTED IN THE MANUAL.

GENERAL SAFETY INSTRUCTIONS & WARNINGS

- If the surface of the scale becomes dirty, wipe it with a clean cloth and alcohol or glass cleaner. Keep it away from soap, water, abrasive liquids, high heat, or extreme cold.
- Check the battery if the scale is not functioning properly, and replace it if necessary.
- Store the batteries in a place out of reach of young children and take precautions to prevent accidental swallowing. Seek medical help immediately if a battery is swallowed.
- If the scale will not be used for a long period, remove the battery and store it separately.
- If the batteries are damaged or leaking, dispose of them safely and avoid contact.

FEATURES

- Maximum weight capacity: 3 kg - 180 kg
- High-precision measurement sensors (Division: 1D = 50 g)
- Unit conversion option (kg/lb/st)
- Auto power-off function
- Low battery and overload indicators

INSTALLING THE BATTERIES

- Open the battery compartment cover.
- Insert 2 AAA 1.5V batteries, ensuring correct polarity, and close the cover.
- When the battery level is low, a warning symbol will appear on the display. Replace the batteries immediately.

INSTRUCTIONS FOR USE

- Place the scale on a flat surface. (Measurements taken on uneven surfaces or soft surfaces like carpets and rugs may be inaccurate.)
- Press the small button located at the bottom of the scale to select the weight unit (kg, lb, or st).
- Step onto the scale barefoot and balance your weight.
- The measurement will appear on the display approximately 3 seconds after stepping on the scale.
- The scale will automatically turn off after 8 seconds.

Notes

- Ensure the scale resets to zero before each use.
- Use the scale only for personal weight measurement.
- Remove the batteries during long-term storage.

MAINTENANCE AND CLEANING

- Wipe the scale with a damp cloth. Do not use detergents or abrasive cleaners.
- Store the device in an upright position to protect the sensitive sensors.

OKOK SETUP

S1 - Install App

- Switch on Bluetooth
 - Phone or Tablet connects to Net (Wi-Fi, 3G or 4G)
- Use your device to scan QR Code or search OKOK International from App Store/Google Play to install App.


S2 - Pair the Devices, Register and Complete User Information

Open APP, begin to weigh and pair.

Register as a user.

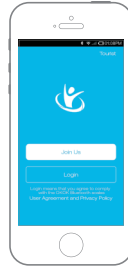Complete personal information.

S3 - Weighing

Enter into "Dynamic Weight" interface.

Stand on the scale and synchronize weight data onto APP.


S4 - View the Test Result

Enter into "Weight Trend" interface.

View test report.

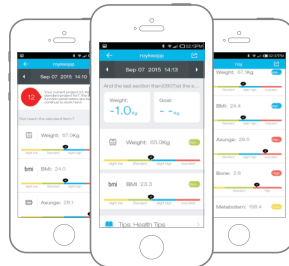Correct Disposal of this product

- This marking indicates that this product should not be disposed with other household wastes throughout the EU.
- To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.
- To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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